

# CIBS OFFICE

## Wireless telephony at work for mobile- and desk telephones


Many people who are hard of hearing find it difficult to keep up with a demanding routine filled with telephones, conversations with clients and colleagues and the use of audio equipment. CIBS OFFICE consists of a multifunctional Tango Bluetooth amplifier and a Maestro Communicator neck loop. It provides a wireless solution that covers the various needs that a person who is hard of hearing may have at work.

CIBS OFFICE gives you the possibility of wireless amplification of mobile-, cordless- and desk phones. Digital signal processing and the possibility of high-quality sound for both ears means that speech recognition and easier handling of telephone conversations again become a reality.

The advantage of Tango Bluetooth is that you can still use the same telephone as your colleagues and keep all the functions you have by using system telephony at work. Most of the desk telephones do not have Bluetooth, so by connecting to Tango Bluetooth you get sound transmitted wirelessly through the Maestro Communicator. If you use mobile phones or wireless phones with Bluetooth, you will receive the telephone sound directly through Maestro.

By purchasing a handset lifter (accessory), it is possible to be notified about and to carry out a phone conversation up to 50-70 meters from the desk phone.

- Digital signal handling provides you with the latest technology within speech recognition and flexibility
- Separate volume and tone adjustment, with LED light that remembers the last setting on each function
- Simple connection and calibration of Tango Bluetooth
- Fits most types of telephones
- Remarkable range: up to 70 meters from Tango Bluetooth

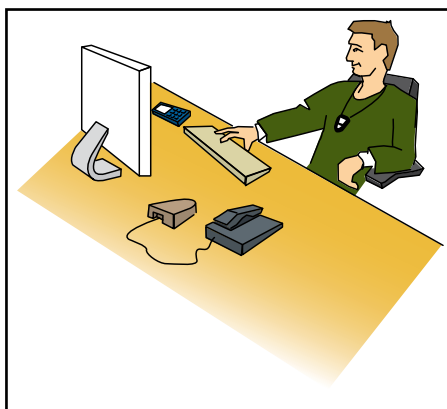

 Bluetooth®

**Wireless solution for your mobile and desk phone**

- Very easy to use
- Integrated speech and audio amplifier
- Active noise reduction
- Stepless volume control
- Rechargeable battery on Maestro Communicator
- Security against tapping
- Simple connection
- No Channels (Like FM)

Art.No: 108-111

**COMMidt®**

Conversations with colleagues, clients and at meeting are often challenging in noisy and busy areas. The Maestro Communicator has a built-in speech amplifier which transfers the wireless digital sound from your dialogue partners. For larger meetings and distances, it may be a good idea to use a Bluetooth microphone such as Allegro or Liberto (see accessories) to ensure that the distance between the sound source and the noisy surrounding is not a problem.

CIBS OFFICE amplifies audio equipment e.g. PC, radio, CD/MP3 players, so that you can enjoy wireless music or PC-sound with crystal clear sound without disturbing others.

- Gives very good and high sound quality. Suitable for all types of hearing aid(s) / CI and hearing losses
- Very suitable for individual adapting and adjusting of workplace
- Choice of Telephone, Neck Loop/Headset and Dialogue functions
- Separate volume and tone adjustments for all functions
- Integrated voice amplifier (Dialogue)
- Adjustable microphone sensitivity
- Simple configuration for most telephones at home and at work
- Stylish design, flexible and easy to use!

#### Inclusive:

Maestro, Tango VII, USB cable, clips, connection cable (cable between Tango VII and the telephone) and international charger

#### Accessories:

Product #	Product Name	Product Description
108-041	Audiostreamer	Provides sound via Bluetooth from TV/HiFi and alerting system to Maestro
108-021	Allegro	Bluetooth microphone with push-to-talk button
108-031	Liberto	Bluetooth combined microphone with omni- and directional zoom functionality and with push-to-talk button
108-051	Mic. Controller	Controls up to two Liberto or Allegro mics. via Bluetooth
227-001	Tango Bluetooth	Amplifier with Bluetooth for desk telephones to be used with Maestro
957-001	Twin charger cradle	Double charger
901-001	External Clip-On Microphone 40 cm	2.5mm jack. Enables you to wear Maestro under clothing
902-001	External Clip-On Microphone 150 cm	2.5mm jack. For conversation use, e.g. in the car
907-001	12/24V Car Charger	To charge Maestro in the car
908-001	Audio Cable for TV/CD/MP3	Enables to listen to audio units without Bluetooth
915-001	Omni Directional Table Microphone	Corded microphone with battery. Inclusive adapter cable
914-002	Directional Microphone	Corded directional microphone
926-001	Scart/Phono Audio Cable	500 cm cable
530-001	Earbuds (Headset) with Microphone	With 3.5 mm jack
521-002	Headset with Microphone	With 3.5 mm jack
527-002	Headset (listen only)	With 3.5 mm jack
930-001	USB Cable	Charge and programme Maestro on the computer
955-002	International Charger (mains power)	
910-001	Handset Lifter	GN Netcom
521-002	Headset with Microphone	With 3.5 mm jack
920-001	Audio Cable, 3,5mm Jack	Connect audio units like MP3, PC, radio, CD player etc
909-001	Audio Extension Cable, 300cm	Extension for the Audio Cable product # 920

#### Tech. info, Maestro

- Bluetooth 2.0
- BlueCore 5 Multimedia Engine
- Talk time: up to 15 hours
- Standby time: > 500 hours
- Charging time: 1.5 hours
- Rechargeable battery
- Battery type: Li-Ion Polymer 3.7 volt
- Dimension: 64 (L) x 35 (W) x 17 (H) mm (2,5 (L) x 1,4 (W) x 0,7 (H) inch.)
- Weight: approx. 40 grams (approx. 1,4 ounces)
- Headset jack (standard 3.5mm Jack)
- Simple unit pairing
- Headset, hands-free and A2DP (Hi-Fi) profile

#### Tech. info, Tango VII

- Bluetooth 2.0
- Class 1 device
- 1 Audio output
- 1 microphone input
- Dimension: 110 (L) x 70 (W) x 50 (H) mm (43 (L) x 27 (W) x 20 (H) inch.)
- Simple unit pairing