

MIC. CONTROLLER


Select the sound source you want to listen to!

Mic.Controller is a remote control for the Bluetooth microphones Liberto and Allegro to be used with the Maestro Communicator. The unit is connected to Maestro via a cable and provides you with the option to choose which microphone(s) you wish to listen to. This provides great flexibility, especially for challenging situations with background noise, like at meetings or social gatherings. The Mic.Controller is easy to use and by using the two control buttons you can decide which microphone you wish to receive sound from.

Easy control provides the options to:

- Listen to one microphone
- Listen to two microphones simultaneously
- Select actively between 2 microphones
- Quiet mode

International charger and connection cable to Maestro are included

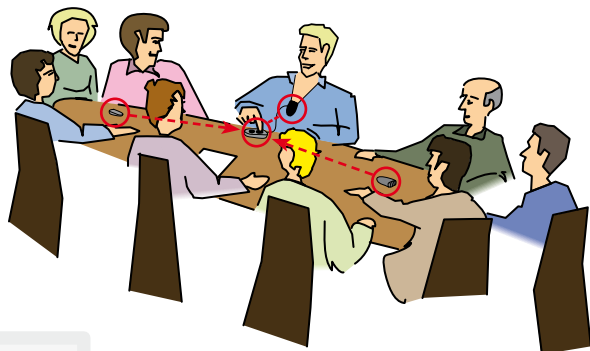

Tech. info

- Bluetooth 2.0
- BlueCore05 Multimedia Engine
- Talk time: up to 15 hours
- Standby time: > 500 hours
- Charging time: 1.5 hours
- Rechargeable battery
- Battery type: Li-Ion Polymer 3.7 volt
- Dimension: 64 (L) x 35 (W) x 17 (H) mm (2,5 (L) x 1,4 (W) x 0,7 (H) inch.)
- Weight: approx. 40 grams (approx. 1,4 ounces)
- Simple unit pairing
- Headset, hands-free and A2DP (Hi-Fi) profile



Scale 1:1

Easy remote control of your microphones

- Freedom to choose the sound source
- Easy control and connection
- Long range: up to 30 metres
- Unique battery capacity
- Very easy to use
- Eavesdrop secure
- Simple connection
- No channels (like FM)

Art.No: 108-051